



sam_2026-10-08_09-37-12_CFX96-HSP90-cGAS.pcrd

10/8/2026 17:36

Report Information

User: BioRad/sam
Data File Name: sam_2026-10-08_09-37-12_CFX96-HSP90-cGAS.pcrd
Data File Path: C:\Users\Samb\Desktop
Well Group Name: All Wells
Report Differs from Last Save: Yes

Run Setup

Run Information

Run Date: 10/8/2026 09:37
Run User: sam
Run Type: User-defined
Plate File: 20261008-mgig-HSP90-cGAS-polyIC-valentina.pltd
ID:
Notes:
Sample Volume: 20
Temperature Control Mode: Calculated
Lid Temperature: 105
Base Serial Number: CC009827
Optical Head Serial Number: 785BR3659

Protocol

- 1: 95.0°C for 0:30
- 2: 95.0°C for 0:03
- 3: 60.0°C for 0:05
- Plate Read
- 4: GOTO 2, 39 more times
- 5: Melt Curve 65.0°C to 95.0°C: Increment 0.5°C 0:05
- Plate Read

Plate Display


	1	2	3	4	5	6	7	8	9	10	11	12
A	Unk-1 HSP90 A1PC	Unk-1 HSP90 A1PC	Unk-1 HSP90 A1PC	Unk-2 HSP90 A2PC	Unk-2 HSP90 A2PC	Unk-2 HSP90 A2PC	Unk-3 HSP90 A3PC	Unk-3 HSP90 A3PC	Unk-3 HSP90 A3PC	Unk-4 HSP90 A4PC	Unk-4 HSP90 A4PC	Unk-4 HSP90 A4PC
B	Unk-5 HSP90 A5PC	Unk-5 HSP90 A5PC	Unk-5 HSP90 A5PC	Unk-6 HSP90 B1PC	Unk-6 HSP90 B1PC	Unk-6 HSP90 B1PC	Unk-7 HSP90 B2PC	Unk-7 HSP90 B2PC	Unk-7 HSP90 B2PC	Unk-8 HSP90 B3PC	Unk-8 HSP90 B3PC	Unk-8 HSP90 B3PC
C	Unk-9 HSP90 B4PC	Unk-9 HSP90 B4PC	Unk-9 HSP90 B4PC	Unk-10 HSP90 B5PC	Unk-10 HSP90 B5PC	Unk-10 HSP90 B5PC	NTC-1 HSP90	NTC-1 HSP90	NTC-1 HSP90			
D	*Unk-11 cGAS A1PC	*Unk-11 cGAS A1PC	*Unk-11 cGAS A1PC	*Unk-12 cGAS A2PC	*Unk-12 cGAS A2PC	*Unk-12 cGAS A2PC	*Unk-13 cGAS A3PC	*Unk-13 cGAS A3PC	*Unk-13 cGAS A3PC	*Unk-14 cGAS A4PC	*Unk-14 cGAS A4PC	*Unk-14 cGAS A4PC

Plate Display

	1	2	3	4	5	6	7	8	9	10	11	12
E	*Unk-15 cGAS A5PC	*Unk-15 cGAS A5PC	*Unk-15 cGAS A5PC	*Unk-16 cGAS B1PC	*Unk-16 cGAS B1PC	*Unk-16 cGAS B1PC	*Unk-17 cGAS B2PC	*Unk-17 cGAS B2PC	*Unk-17 cGAS B2PC	*Unk-18 cGAS B3PC	*Unk-18 cGAS B3PC	*Unk-18 cGAS B3PC
F	*Unk-19 cGAS B4PC	*Unk-19 cGAS B4PC	*Unk-19 cGAS B4PC	*Unk-20 cGAS B5PC	*Unk-20 cGAS B5PC	*Unk-20 cGAS B5PC	*NTC-2 cGAS	*NTC-2 cGAS	*NTC-2 cGAS			
G												
H												

Quantification

Step #: 3
Analysis Mode: Target
Cq Determination: Single Threshold
Baseline Method:
HSP90: Auto Calculated
Threshold Setting:
HSP90: 4696.48, Auto Calculated

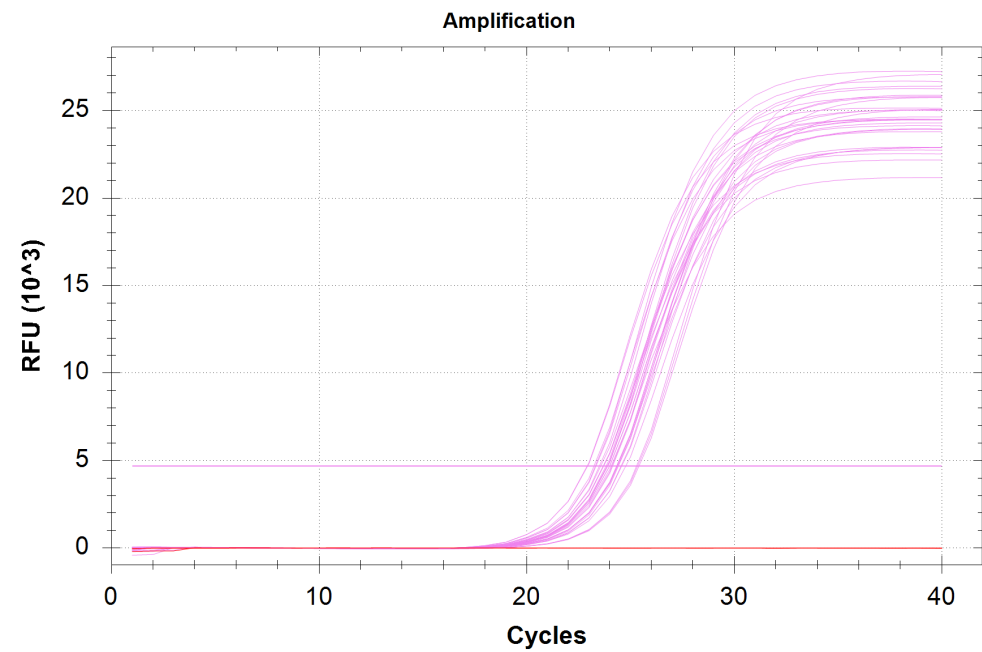

Quantification Data

Well	Fluor	Target	Content	Sample	Cq	Cq Mean	Cq Std. Dev
A01	SYBR	HSP90	Unkn-01	A1PC	24.39	24.36	0.028
A02	SYBR	HSP90	Unkn-01	A1PC	24.36	24.36	0.028
A03	SYBR	HSP90	Unkn-01	A1PC	24.33	24.36	0.028
A04	SYBR	HSP90	Unkn-02	A2PC	23.88	23.76	0.099
A05	SYBR	HSP90	Unkn-02	A2PC	23.71	23.76	0.099
A06	SYBR	HSP90	Unkn-02	A2PC	23.70	23.76	0.099
A07	SYBR	HSP90	Unkn-03	A3PC	23.89	23.97	0.063
A08	SYBR	HSP90	Unkn-03	A3PC	24.00	23.97	0.063
A09	SYBR	HSP90	Unkn-03	A3PC	24.00	23.97	0.063
A10	SYBR	HSP90	Unkn-04	A4PC	23.57	23.78	0.264

Quantification Data

Well	Fluor	Target	Content	Sample	Cq	Cq Mean	Cq Std. Dev
A11	SYBR	HSP90	Unkn-04	A4PC	23.69	23.78	0.264
A12	SYBR	HSP90	Unkn-04	A4PC	24.08	23.78	0.264
B01	SYBR	HSP90	Unkn-05	A5PC	25.41	25.35	0.059
B02	SYBR	HSP90	Unkn-05	A5PC	25.35	25.35	0.059
B03	SYBR	HSP90	Unkn-05	A5PC	25.29	25.35	0.059
B04	SYBR	HSP90	Unkn-06	B1PC	23.51	23.39	0.105
B05	SYBR	HSP90	Unkn-06	B1PC	23.31	23.39	0.105
B06	SYBR	HSP90	Unkn-06	B1PC	23.34	23.39	0.105
B07	SYBR	HSP90	Unkn-07	B2PC	24.55	24.44	0.104
B08	SYBR	HSP90	Unkn-07	B2PC	24.42	24.44	0.104
B09	SYBR	HSP90	Unkn-07	B2PC	24.35	24.44	0.104
B10	SYBR	HSP90	Unkn-08	B3PC	24.58	24.50	0.323
B11	SYBR	HSP90	Unkn-08	B3PC	24.14	24.50	0.323
B12	SYBR	HSP90	Unkn-08	B3PC	24.78	24.50	0.323
C01	SYBR	HSP90	Unkn-09	B4PC	23.25	23.03	0.183
C02	SYBR	HSP90	Unkn-09	B4PC	22.92	23.03	0.183
C03	SYBR	HSP90	Unkn-09	B4PC	22.94	23.03	0.183
C04	SYBR	HSP90	Unkn-10	B5PC	23.82	23.95	0.149
C05	SYBR	HSP90	Unkn-10	B5PC	23.91	23.95	0.149
C06	SYBR	HSP90	Unkn-10	B5PC	24.11	23.95	0.149
C07	SYBR	HSP90	NTC-01		N/A	0.00	0.000
C08	SYBR	HSP90	NTC-01		N/A	0.00	0.000
C09	SYBR	HSP90	NTC-01		N/A	0.00	0.000

QC Parameters

Data

Description	Value	Use	Results	Exclude Wells	All excluded wells
Negative control with a Cq less than	38	True		False	
NTC with a Cq less than	38	True		False	
NRT with a Cq less than	38	True		False	
Positive control with a Cq greater than	30	True		False	
Unknown without a Cq	N/A	True		False	
Standard without a Cq	N/A	True		False	

Data

Description	Value	Use	Results	Exclude Wells	All excluded wells
Efficiency greater than	110.0	True			
Efficiency less than	90.0	True			
Std Curve R^2 less than	0.980	True			
Replicate group Cq Std Dev greater than	0.20	True	HSP90:A10, A11, A12, B10, B11, B12.	False	